

ORIGINAL RESEARCH

SERUM β_2 -MICROGLOBULIN VALUES AMONG HEALTHY BRAZILIANS USING A DPC IMMULITE® ASSAY

Fabiola Branco Filippin and Liliete Canes Souza

FILIPPIN FB et al. Serum β_2 -Microglobulin values among healthy Brazilians using DPC IMMULITE® Assay. **CLINICS** 60(1):47-50, 2005.

PURPOSE: The present study was designed to determine the normal range of serum beta 2-microglobulin (Sb₂M) levels among healthy volunteers in Brazil. Levels of Sb₂M are elevated in human immunodeficiency virus (HIV)-infected patients and have been shown to be the best predictor of HIV infection status and of some malignant disorders, especially multiple myeloma. In order to achieve its optimal use in Brazilian clinical diagnosis, an adequate reference interval study was performed for Sb₂M IMMULITE® assay, based on the fact that its reference range limits were evaluated among European populations.

METHODS: Ninety-six healthy blood donors were evaluated, and Sb₂M levels were measured by chemiluminescent enzyme immunoassay using the IMMULITE® automated analyzer.

RESULTS: A normal range of Sb₂M values, established by a nonparametric statistical method, was 1.05 to 3.9 mg/mL, with the upper limit being higher than that reported elsewhere.

CONCLUSIONS: This study presented new data indicating that there is a significant difference between the current reference limits for Sb₂M IMMULITE® assay and those found in Brazil, providing evidence that significant differences in range of normal values may occur among different populations, and that these new values should be considered for Brazilian individuals.

KEYWORDS: β_2 -microglobulin. Reference interval. Chemiluminescent immunoassay. Healthy volunteers. Prognostic value.

Beta 2-microglobulin (b2M) a protein found in all nucleated cells, is associated with histocompatibility class 1 antigens on cell surface membranes (particularly abundant on lymphocytes and monocytes).¹ It was first discovered in the urine of patients with renal failure.² Because of its low molecular weight (11,800 daltons), 95% of all free b2M in plasma is eliminated by glomerular filtration, and 99.9% of it taken up by proximal tubular cells. In the presence of a normal renal filtration rate, elevated serum b2M (Sb₂M) levels indicate high b2M production or release.³ Upon ac-

tivation of the immune system, both B- and T-lymphocytes actively release b2M into circulation from where it is later eliminated via glomerular filtration and tubular reabsorption.⁴ Increased Sb₂M concentrations have been found in patients with blood cells dyscrasias, including multiple myeloma,^{5,6} malignant lymphoproliferative disorders,⁷ myeloproliferative disorders,⁸ certain viral infections such as HIV infection⁹, cytomegalovirus, non-A and non-B hepatitis, and infectious mononucleosis.¹⁰ The level of Sb₂M is also used as an important prognostic factor in patients with multiple myeloma^{5,6} and HIV infections.⁹

The normal reference range for Sb₂M in healthy individuals found in the literature was determined using North American and European populations (DPC Immulite® assay). In an attempt to determine the normal range of Sb₂M levels among

From the Department of Clinical Analysis, Science and Health Center, Federal University of Santa Catarina - Florianópolis/SC, Brazil.

E-mail: liliete@ccs.ufsc.br

Received for publication on September 01, 2004.

Accepted for publication on October 15, 2004.

This is an Open Access article distributed under the terms of the Creative Commons Attribution Non-Commercial License (<http://creativecommons.org/licenses/by-nc/3.0/>) which permits unrestricted non-commercial use, distribution, and reproduction in any medium, provided the original work is properly cited.

asymptomatic seronegative African individuals, Piwowar et al. (1995)¹⁰ found that 50% of these subjects had abnormal values (mean values 2.35 mg/mL). In another study, the mean Sb2M value among healthy volunteers was 1.36 mg/mL.¹¹ In this report, we describe Sb2M values in healthy Brazilian subjects from a group of blood donors at the Hematology and Hemotherapy Center of Santa Catarina (*HEMOSC*).

SUBJECTS

Blood donors, men and nonpregnant women ranging from 25 to 50 years of age, were included in the study. All individuals underwent a clinical assessment evaluated by *HEMOSC*, including questionnaires and screening laboratory tests according to Brazilian legal criteria. All volunteers fulfilling these criteria gave verbal and written consent to have an additional 5 mL of blood drawn for the study purpose. The experimental protocol was approved by the Ethical Committee of the University Hospital of Florianópolis, Federal University of Santa Catarina.

METHODS

Whole blood was collected with a Vacutainer system in 10-mL dry tubes, and serum was separated by centrifugation. In order to eliminate lipid components, serum samples were ultracentrifuged at 18000 g for 20 min. All blood samples were obtained between 2:00 p.m. and 6:00 p.m. and were processed within 5 h after venipuncture. Samples were prediluted 1-in-41 in b2M Sample Diluent (IMMULITE® Beta-2 Microglobulin kit, DPC, USA) according to the manufacturer's instructions, and 5 mL were required for the assay procedure. Sb2M was measured by a solid-phase 2-

site chemiluminescent enzyme immunoassay using the IMMULITE® Automated Analyzer (DPC, Los Angeles, USA).

Statistical analysis

In the reference range study for the DPC IMMULITE® assay, the mean value was 1.26 mg/mL. The National Committee for Clinical Laboratory Standards (NCCLS) addresses the use of a priori and a posteriori exclusion criteria to obtain a healthy samples.¹² Therefore, all subjects with Sb₂M values greater than 2.0 mg/mL were clinically reevaluated for value acceptance to ensure that every individual was "healthy". Data were entered and analyzed with Microsoft Excel software. The mean and standard deviation values were calculated. The 95th-percentile reference ranges were determined by using 2.5 and 97.5 percentiles. According to the procedure recommended by the NCCLS, the observations were ranked according to size, and the 2.5 and 97.5 percentiles were obtained as the 0.025 (n + 1) and 0.975 (n + 1) ordered observations respectively (n= sample size). If the estimated rank values were not integers, linear interpolation was carried out.¹²

RESULTS

Ninety-six Brazilian blood donors, 60 men and 36 women, with mean age (yr) of 33, were included in the study. Values of Sb2M ranged between 1.05 and 3.9 mg/mL, as shown in Figure 1, and the mean was 2.46 mg/mL (standard deviation [SD] = 0.73). Therefore, using 97.5th percentile according to a routine nonparametric method,¹² 3.9 mg/mL was considered the upper normal limit. The range established was 1.3–3.9 mg/mL.

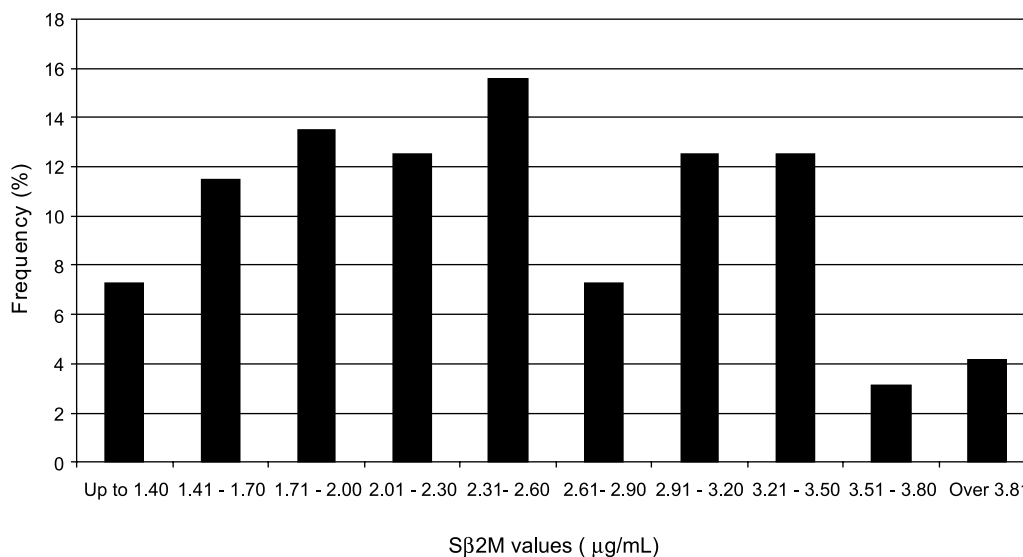


Figure 1 - Serum beta 2-microglobulin (Sb2M) levels in blood donors (mg/mL) plotted versus its frequency distribution in percentage.

DISCUSSION

Researchers worldwide have reported the usefulness of Sb2M as a reliable maker in evaluating progression of HIV-infection and multiple myeloma.^{5-8,10,13} In HIV-infected subjects, Sb2M appears to be a better predictor of HIV infection status than CD4 counts and also a better predictor of survival, with increased Sb2M levels being associated with disease progression.¹⁰

Increased production and consequent release of Sb2M into body fluids have been reported in patients with multiple myeloma. The importance of these data for clinicians is based on the fact that myeloma plasma cells produce Sb2M. Therefore, it is assumed that Sb2M values could indicate not only the size but also the proliferative activity of the myeloma cells¹³ and could be accepted as a nonspecific maker of tumor mass.^{5,6,13,14} In addition, the prognostic role of Sb2M in other proliferative disorders has been confirmed by many authors,^{7,14,15} facilitating clinical assessment.¹⁵ On the other hand, it has been reported that in patients with low glomerular filtration rate, Sb2M levels were raised and laboratory findings should be interpreted with caution.^{6,7,13}

Reference intervals are used clinically, together with additional information, as guidelines concerning the state of the patient. The reference range study for IMMULITE® Sb2M (1.0–1.7 mg/mL, mean 1.26 mg/mL) defined by the manufacturer was determined using serum samples from almost 800 adult volunteers, including both men and non-pregnant women, ranging from 20 to 70 years of age. The subjects were selected based on questionnaires. Samples were collected in France, Germany, the Netherlands, and Portugal, and values obtained were analyzed

nonparametrically. The suggested upper limit was 1.7 mg/mL. These values may not be representative for Brazilian individuals. In the absence of adequate reference range limits for Sb2M derived from any Brazilian population, many investigators have been using data obtained from North American and European subjects that presented upper limits 2.3 times lower than that obtained in the present study.

The NCCLS recommends that nonparametric reference intervals be used and that the sample size, *n*, should consist of at least 120 values in order to compute a 95% reference interval, 2.5% and 97.5% points of the distribution¹². This number is the minimum number of values needed for construction of the confidence intervals giving information about the accuracy of the calculated reference limits.

In this study, 96 samples were evaluated. Our goal was not to produce new reference intervals, but to bring more focus on the need for every laboratory to produce its own reference intervals for all parameters in spite of many difficulties, including practical ones and high expenses.

The present study indicates that significantly different ranges of Sb2M values may occur among different populations, suggesting that each Brazilian laboratory should establish its own reference limits.

ACKNOWLEDGEMENTS

We thank the Hematology and Hemotherapy Center of Santa Catarina and the Clinical Analysis Laboratory, University Hospital of Florianópolis for the samples and technical assistance. We are grateful to DPC IMMULITE® for providing the reagent set free of charge for this study and to Dr. Joanita Del Moral from the Hematology Division, University Hospital of Florianópolis.

RESUMO

FILIPPIN FB e col. Valores de β_2 -Microglobulina em Brasileiros saudáveis usando o sistema automatizado DPC IMMULITE®. **CLINICS** 60(1):47-50, 2005.

OBJETIVO: O objetivo do presente estudo foi determinar valores de referência de β_2 -Microglobulina sérica (Sb₂M) em voluntários saudáveis. Sabe-se que tal parâmetro apresenta-se elevado em pacientes infectados com o vírus da imunodeficiência humana (HIV) e tem se mostrado melhor marcador da infecção por HIV e de desordens malignas, especialmente mieloma múltiplo. De forma a se obter o melhor diagnóstico clínico, um intervalo de referência adequado de Sb₂M para a população brasileira

foi determinado empregando-se o ensaio IMMULITE®; já que este tem como parâmetro uma faixa de referência determinada a partir de populações européias.

MÉTODOS: Noventa e seis doadores de sangue saudáveis foram avaliados e os valores de Sb₂M foram medidos por método enzima imunoensaio quimiluminescente usando analisador automatizado IMMULITE®.

RESULTADOS: Os valores de Sb₂M, estabelecidos por método estatístico não paramétrico, apresentaram-se entre 1.05 e 3.9 mg/ml, sendo que o limite superior obtido foi maior que o relatado em outros trabalhos.

CONCLUSÃO: Este estudo apresentou novos valores,

indicando que existe uma diferença significativa entre os limites de referência de Sb2M disponibilizados para o IMMULITE® e os encontrados no presente trabalho, evidenciando a ocorrência de variações importantes entre diferentes populações e que novos valores devem ser

considerados para brasileiros.

UNITERMOS: β_2 -microglobulina. Intervalo de referência. Imunoensaio quimiluminescente. Voluntários saudáveis. Valor prognóstico.

REFERENCES

1. Peterson PA, Cunningham BA, Edelman GM. 2-Microglobulin: a free immunoglobulin domain. *Proc Natl Acad Sci USA* 1972;69:1697-701.
2. Berggard I & Bearn AG. Isolation and properties of a low molecular weight beta-2-globulin occurring in human biological fluids. *J Biol Chem* 1968;243(15):4095-103.
3. Schardijn GHC & van Eps LWS. Beta 2-microglobulin: its significance in the evaluation of renal function. *Kidney Int* 1987;32(5):635-41.
4. Schweizer RT, Moore R, Bartus SA, Bow L & Hayden J. Beta 2-microglobulin monitoring after renal transplantation. *Transplant Proc* 1981;13(3):1620-23.
5. Durie BG, Stock-Novak D, Salmon SE, Finley P, Beckord J, Crowley J. Prognostic value of pretreatment serum beta 2 microglobulin in myeloma: a southwest oncology group study. *Blood* 1990;75(4):823-30.
6. Greipp PR, Katzman JA, O'Fallon WM, Kyle RA. Value of beta 2-microglobulin level and plasma cell labeling indices as prognostic factors in patients with newly diagnosed myeloma. *Blood* 1988;72(1):219-23.
7. Constantines IP, Pathouli C, Karvountzis G, Papadopoulos P, Varvoutsi-Constantinides M, Eliakis P, et al. Serum beta 2 microglobulin in malignant lymphoproliferative disorders. *Cancer* 1985;15(10):2384-89.
8. Rodriguez J, Cortes J, Talpaz M, O'Brien S, Smith TL, Rios MB, et al. Serum beta-2 microglobulin levels are a significant prognostic factor in Philadelphia chromosome-positive chronic myelogenous leukemia. *Clin Cancer Res* 2000;6:147-52.
9. Piwowar EM, Tugume SB, Grant RM, Lutalo T, Pattishall K, Katongole-Mbidde E. Beta-2 microglobulin values among human immunodeficiency virus (HIV)-negative, HIV-positive asymptomatic, and HIV-positive symptomatic Ugandans. *Clin Diagn Lab Immunol* 1995;2:236-7.
10. Cooper E, Forbes M & Hamblin M. Serum beta 2-microglobulin and C-reactive protein concentrations in viral infections. *J Clin Pathol* 1985;37:1140-3.
11. Diamondstone LS, Tollerud DJ, Fuchs D, Wachter H, Brown LM, Maloney E, et al. Factors influencing serum neopterin and beta 2-microglobulin levels in a healthy diverse population. *J Clin Immunol* 1994;14(6):368-74.
12. National Committee for Clinical Laboratory Standards. How to define and determine reference intervals in the clinical laboratory: approved guideline, NCCLS document C28-A and C28-A2, Villanova, PA, 1995, 2001, NCCLS.
13. Bataille R, Greinier J, Sany J. Beta-2-microglobulin in myeloma: optimal use for staging, prognostic, and treatment—a prospective study of 160 patients. *Blood* 1984;63(2):468-76.
14. Mellito L, Musto P, Tomasi P, Cascavilla N, Bodenizza C, Ladogana S, et al. Serum beta 2-microglobulin in malignant lymphoproliferative disorders. *Tumori* 1988;74(2):129-35.
15. Dimopoulos MA, Cabanillas F, Lee JJ, Swan F, Fuller L, Allen PK, et al. Prognostic role of serum beta 2-microglobulin in Hodgkin's disease. *J Clin Oncol* 1993;11(6):1108-11.