



Cytotoxicity assessment of three endodontic sealing cements used in periapical surgery. *In vitro* study

Evaluación de la citotoxicidad de tres cementos selladores endodóncicos utilizados en cirugía periapical: estudio in vitro

Mildred Martínez-Cortés,* Carlos Tinajero-Morales,* Carlos Rosales,§ Eileen Uribe-Querol^{||}

ABSTRACT

Introduction: Presently there are many retrofilling materials in the market, nevertheless, little is known about their toxicity on gingival fibroblasts. **Objective:** To assess cytotoxicity of three materials to human gingival fibroblasts and L929 mouse fibroblasts cell line. **Material and methods:** EndoSequence[®] BC RRM[™] (ERRM; root repair material), white MTA Angelus[®] (MTA) and intermediate restoration material (IRM[®]) conditioned media were obtained when materials were freshly mixed, at setting time and after 1, 24 and 72 hours of setting time. Cell morphology was assessed with light microscopy and cell viability was assessed through mitochondrial metabolic activity with 3-(4,5-dimethylthiazole-2-yl)-2,5-diphenyl tetrazolium bromide (MTT). Statistical analysis was conducted with ANOVA. **Results:** We found that ERRM material did not exhibit cytotoxic effects on used fibroblast, nevertheless, MTA and IRM[®] respectively exhibited moderate and severe cytotoxicity, thus indicating the materials were not fully harmless. **Conclusion:** Bioceramic cements like ERRM could be considered the most compatible retrofilling-materials.

Key words: Cell morphology, cytotoxicity, retro-filling, endodontic sealers, MTT, inflammation.

Palabras clave: Morfología celular, citotoxicidad, retroobturación, selladores endodóncicos, MTT, inflamación.

RESUMEN

Introducción: Existen diversos materiales de retroobturación, pero poco se sabe de su toxicidad sobre fibroblastos gingivales. **Objetivo:** Evaluar la citotoxicidad de tres materiales de retroobturación sobre fibroblastos gingivales humanos y fibroblastos de la línea L929. **Material y métodos:** Los medios condicionados de los materiales de retroobturación EndoSequence[®] BC RRM[™] (ERRM), trióxido mineral agregado MTA Angelus[®] blanco (MTA) y material de restauración intermedia (IRM[®]) se obtuvieron en fresco, al tiempo de fraguado, y después de 1, 24 y 72 horas del tiempo de fraguado. La morfología celular fue evaluada por microscopía de luz y la viabilidad celular fue evaluada a través de la actividad metabólica mitocondrial con 3-(4,5-dimetiltiazol-2-il)-2,5-difenil bromuro de tetrazolio (MTT). El análisis estadístico se realizó por ANOVA. **Resultados:** El material ERRM no mostró efectos citotóxicos sobre los fibroblastos. Sin embargo, el MTA y el IRM[®] mostraron citotoxicidad moderada y alta, respectivamente. Esto revela que el MTA y el IRM[®] no son completamente inocuos. **Conclusión:** Los materiales biocerámicos como el ERRM pueden ser considerados los materiales de retroobturación más biocompatibles.

INTRODUCTION

In the field of Endodontics, surgery is recommended in cases when canal treatment failed or when there is any contraindication to conducting conventional root canal treatment.¹ In an event of periapical surgery, an incision is first performed in order to uncover the apical section of the infected tooth. After this, infected tissue in the affected area is removed with a curette. Retrofilling site is later prepared, retrofilling is performed and, finally, incision is sutured. In this procedure, retrofilling materials are used which ideally would provide apical seal and facilitate periapical tissue reparation.² Microscope use, surgical techniques developments, improvement of lighting and availability of new retrofilling materials greatly contribute to the success of periapical

* Department of Endodontics, División de Estudios de Postgrado e Investigación.

§ Immunology Department, Instituto de Investigaciones Biomédicas.

^{||} Neurobiology and Development, División de Estudios de Postgrado e Investigación.

Universidad Nacional Autónoma de México, Mexico City, Mexico.

Received: January 2016. Accepted: October 2016.

© 2016 Universidad Nacional Autónoma de México, [Facultad de Odontología]. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

This article can be read in its full version in the following page: <http://www.medigraphic.com/facultadodontologiaunam>

surgery.^{3,4} These advantages have heralded a new era in endodontic surgery.⁵

An ideal retrofilling material must possess suitable adhesion to dentin, appropriate sealing, dimensional stability, insolubility to tissue fluids, fair compressibility, adequate working time and rapid setting. Moreover, it must be resorbable, radio-opaque, easy to manipulate and biocompatible with the host's tissues.⁶⁻⁸ Among materials used in retrofilling we can count the following: amalgam, resin materials, zinc oxide-eugenol cements, glass ionomer cements, polycarboxylate cements Cavit™ and gutta-percha.⁹

Amalgam was the most frequently used material in periapical surgeries, nevertheless, it did not provide suitable clinical results.⁴ After amalgam, use of cements based on zinc oxide and eugenol was established. Some of the disadvantages of these cements are long setting time and high solubility. Incorporation of polymethylmethacrylate was an improvement in these cements. Intermediate restoration material (IRM®) is the outcome of this improvement. This cement is easy to handle, easy to mix, possesses better sealing capacity, requires lesser setting time⁴ and provides suitable post-operative results.¹⁰ Nevertheless, this cement is not biocompatible, since it elicits inflammatory response in tissues surrounding the tooth.³ Another cement used is mineral trioxide aggregate (MTA), composed mainly of calcium silicate. It was developed at Loma Linda University, California, in 1993, later available in the market from 1999 onwards. MTA possesses excellent sealing capacity and provides suitable long term results. Moreover, it is a tissue-biocompatible material since it does not cause inflammatory reactions around teeth.¹¹ Furthermore, it possesses reparation capacity through hard tissue formation and counts with antibacterial properties.⁷ Although MTA as retrofilling material has won recognition from the time it was introduced in the endodontic market, it is a material difficult to handle¹² and hard to place during treatment.¹³ An additional drawback of this material is its prolonged setting time.^{6,11} MTA has been improved with variations in its composition. MTA Angelus® (Angelus, Londrina Brazil) is one of these variations.¹² White MTA Angelus® has shown to possess greater amounts of calcium carbonate, calcium silicate, and barium and zinc phosphates. These elements contribute to improve setting time and manipulation ease.^{6,14} MTA Angelus® is similar to original MTA with respect to biocompatibility and antimicrobial capacity as well as its physical and chemical characteristics.¹⁵ A new retrofilling material called EndoSequence® BM RRM™ (ERRM, root repair material) is composed of a

combination of calcium silicate and calcium phosphate, and provides biomedical and dental use.¹⁶ It exhibits the following physical properties: it is hydrophilic and homogeneous, since it is presented in a pre-mixed paste, and thus possesses exceptional dimensional stability.^{17,18} Additionally, it is a highly radio-opaque material. This characteristic makes it easy to place during treatment and later to be identified in X-rays. Its 12.8 pH is partly responsible for its anti-bacterial nature, which constantly decreases during a seven day period rendering it thus highly biocompatible.^{6,19} The manufacturer of ERRM claims that it possesses physical and mechanical properties similar to MTA but with better fit and handling characteristics.¹

Biocompatibility

MTA is a biocompatible material.^{8,14} Torabinejad and his team showed that MTA, both freshly mixed and fully set, is less cytotoxic than IRM®.^{8,20} Moreover, MTA does not interfere with cell adhesion.^{8,17} MTA and IRM® have been widely researched and are used as retrofilling materials.¹⁰ Prospective clinical studies have shown that both possess suitable ability for apical healing.³ The manufacturer of ERRM claims that this material possesses excellent physical characteristics as well as suitable biocompatibility, providing very favorable results when it is used as retrofilling material, although there is scant literature supporting this claim. Different retrofilling materials exist in the market, but few comparative reports about their cytotoxicity potential using fibroblasts are available.¹³ Due to the aforementioned facts the aim of the present research paper was to compare cytotoxicity of White MTA Angelus®, IRM® and recently launched in the Mexican market, EndoSequence® BC RRM™ (ERRM) retrofilling materials on human gingival fibroblast culture (FGH) and L929 mouse fibroblast cell line.

MATERIAL AND METHODS

Retrofilling materials

Three retrofilling materials were used: intermediate restoration material (IRM®) which is a cement based on zinc oxide and eugenol, white MTA Angelus® (MTA) and EndoSequence® BC RRM™ root repair material (ERRM) which is a bioceramic material (*Tables I and II*). Materials were prepared as presented in *table II*.

Conditioned media of retrofilling materials

Conditioned media and their respective dilutions were obtained following the norm ISO 10993-5:2009

Table I. Retrofilling material, composition and manufacturer.

Material	Composition	Manufacturer
IRM®	Powder: zinc oxide and polymethylmetacrylate Liquid: eugenol	Dentsply International, Milford DE, USA (Intermediate Restorative Material; Dentsply Sirona Global Headquarters, York, USA)
White MTA Angelus®	Tricalcium silicate, dicalcium silicate, bismuth oxide. Insoluble residue of crystalline silica, calcium oxide and sodium and potassium sulphates	Angelus, Londrina-PR-Brasil
EndoSequence® BC RRM™ (ERRM)	Calcium silicate, zirconium oxide, tantalum pentoxide, calcium monobasic phosphate and filling materials	Brasseler, Savannah GA, USA

Table II. Material preparation.

Sealing material	Preparation	Setting time
IRM®	1.5 spoonful of powder and 1 spoonful liquid	10 minutes
White MTA Angelus®	One spoonful powder one drop distilled water	2.45 hours
EndoSequence® root repair material	Premixed material	2 hours

(International Standardization Organization 2009).²⁰ Briefly, sterile teflon discs were placed on top of sterile Mylar Teflon sheets squares (1.1 x 1.1 cm; both from Tetraflon de México S.A. de C.V. Toluca, Estado de México). Teflon discs were filled with fresh materials to be later placed in tissue culture 24-well plates (Corning Inc., Corning, NY). To this mix, 1.5 mL of supplemented Dulbecco's modified Eagle medium were added (DMEM; GIBCO, Invitrogen, Grand Island, New York, USA). Media was supplemented with 10% bovine fetal serum. (ByProductos SA de CV, Guadalajara, Jalisco, México), 8 mM L-glutamine and antibiotics (10,000 units/mL penicillin and 25 µg/mL Fungizone® (GIBCO, Invitrogen, Grand Island, NY, USA). Culture medium was left in contact with retrofilling materials at different times: fresh (0 h), setting time (*Table II*) and at 1, 24 and 72 hours after setting. Conditioned media was collected in Eppendorf tubes and stored at -20 °C, until used. Media, which were not exposed to any cement, were also collected as control conditioned media, during the same times. Three dilutions of each conditioned medium (1/10, 1/100 and 1/1,000) were prepared.

Cell lines

L929 mouse fibroblast cell line, as well as human gingival fibroblasts (HGF) were used in the present study. HGF were obtained from gingival tissue donated by healthy patients. These procedures were approved by the corresponding Ethics Committee; all donors signed informed consent forms. Gingival

fibroblasts were harvested as follows: gingival tissue was placed in DMEM under sterile conditions, human gingival tissue was washed and divided into small pieces (1 mm³). Gingival tissue pieces were placed in 3 mL DMEM in 50 mL tubes. After this, the medium was removed and new medium with 5 mL of 5 mg/mL collagenase type IV was placed (Sigman-Aldrich Co., St Louis, MO, USA). The tube was placed in water at 37 °C for 2 hours, and it was shaken for 30 seconds every 30 minutes. After this, the tube was centrifuged at 377 g for 2 minutes, supernatant was then removed, and cells were re-suspended in 5 mL supplemented; the medium was then transferred to a 25 cm² culture flask (Corning Inc., Corning NY, USA). The flask was incubated at 37 °C with 5% CO₂ for 24 hours. Non-adhered cells were removed and human gingival fibroblasts were left to grow until reaching confluence.²¹

Cell culture

Both fibroblast cell lines were cultured in supplemented DMEM. Cell cultures were kept in 75 cm² cell culture flasks (Corning Inc., Corning NY, USA) until reaching 80-90% confluence, at which point they were sub-cultured. To this effect, fibroblasts were detached with 0.05% trypsin/EDTA (GIBCO, Invitrogen, Grand Island, NY). Trypsin was neutralized with an equal volume of fresh culture medium. Cell suspension was centrifuged for three minutes at 1,500 rpm. Supernatant was discarded and cell pellet was re-suspended in 5 mL of fresh

medium. One fourth of the re-suspension was placed in a new culture flask. All experiments were conducted between subcultures three and 10.

Cell morphology

In order to assess cell morphology, 10^4 fibroblasts/well were seeded in 96-well plates (Corning Inc., Corning, NY) in a final volume of 100 μ L culture medium per well. Cultures were kept in an incubator for 24 h (Nuairé TM Plymouth, USA) at 37 °C, 95% relative humidity and 5% CO₂. After 24 h, the time required to allow fibroblasts to adhere to the wells, culture medium was removed and 100 μ L of all corresponding conditioned media were placed for another 24 h. Morphology was then evaluated with a model IX70 Olympus microscope (Center Valley, PA, USA). Images were captured with an Evolution-VF Cooled Color camera (Media Cybernetics, Rockville MD, USA) and with the computer program Capture pro 6.0 QImaging (Surrey, British Columbia, Canada).

Metabolic activity assay (MTT reduction)

Cytotoxic potential of retrofilling materials was assessed based on norm ISO 10993-5:2009 (International Standardization Organization) using the 3-(dimethylthiazol-2-yl)-2-diphenyl tetrazolium bromide assay (MTT) (Sigma Aldrich St Louis, MO, USA).²¹ Briefly, the same cultures that were used to take fibroblast pictures were removed from their conditioned media and 50 μ L of dissolution of 1 mg/mL MTT were added into each well. Cultures were protected from the light and kept at 37 °C with 5% CO₂ for 2 h. After this time, MTT was removed and 100 μ L of isopropanol were added (TECSIQUIM, TSQ, Iztacalco, Mexico City). Cultures were incubated for 30 min at room temperature. After this time, dissolution absorbance was measured at 590 nm in a plate reader model Synergy HT, Bio-Tek brand (Vermont USA). Absorbance values were normalized considering 100% as the absorbance obtained from untreated cultures.

Statistics

A one-way ANOVA variance analysis for equal samples was used to determine possible differences among cement groups and among different times. Tukey test was used as post-hoc test, using KaleidaGraph® computer program, Mac version 3.6.2 (Synergy Software, Reading, PA, USA). Conditions were considered statistically different when p value was ≤ 0.05 .

RESULTS

Cells treated with Intermediate restoration material (IRM®) were severely damaged at 24 hours after material had set, but recovered their morphology after 72 hours after the material had set.

Exposing human gingival fibroblasts (*Figure 1A*) or L929 fibroblasts (*Figure 1B*) to IRM® did not affect cell viability when the sealer was fresh, at the setting time or after 1 hour after material was set. However, cell viability was reduced to about 70%, when the conditioning media was from 24 hour after material had set. Cell viability recovered to 50% when fibroblasts were exposed to conditioning media from 72 hours after sealer had set. Cell viability remained high 80 to 90% at all dilutions at all times of conditioning medium from all but 24 hours after sealer had set (*Figures 1A and 1B*).

Treated human gingival fibroblasts (*Figure 1C*) or L929 fibroblasts (*Figure 1D*) had similar morphology to untreated fibroblasts, when exposed to IRM®, freshly mixed, just set and 1 hour after setting. However, by 24 hours after setting, fibroblasts had lost their capacity to maintain a homogenous monolayer (*Figures 1C and 1D* panel e and n, respectively). By 72 hours after setting IRM® caused severe damage L929 cell line, cells had rounded morphology (*Figure 1D* panel o). Cell morphology was not affected with IRM® dilutions (*Figures 1C and 1D* panel g-i and p-r, respectively).

Cells treated with White MTA Angelus® (MTA) remained viable and exhibited typical fibroblasts morphology.

Exposing fibroblasts to MTA did not affect cell viability of either human gingival fibroblast (*Figure 2A*) or L929 cells (*Figure 2B*). Moreover, fibroblasts exposed to serial dilutions of the MTA conditioned medium showed an increased in metabolic activity. Accordingly, metabolic activity was higher than untreated fibroblast after 72 hours after setting at 1/100 and 1/1,000 dilutions (*Figures 2A and 2B*). In agreement, human gingival fibroblasts had a typical bipolar elongated shape and remained attached to the culture dish at all times (*Figure 2C*). Similarly, L929 fibroblasts had a typical triangular elongated shape and remained viable (*Figure 2D*).

Cells treated with EndoSequence® BC RRM™ root repair material (ERRM) remained viable, exhibited typical fibroblasts morphology and displayed higher metabolic activity.

Similarly to MTA, exposing fibroblasts to ERRM did not affect cell viability of either human gingival fibroblast (*Figure 3A*) or L929 cells (*Figure 3B*). Interestingly, fibroblasts exposed to serial dilutions

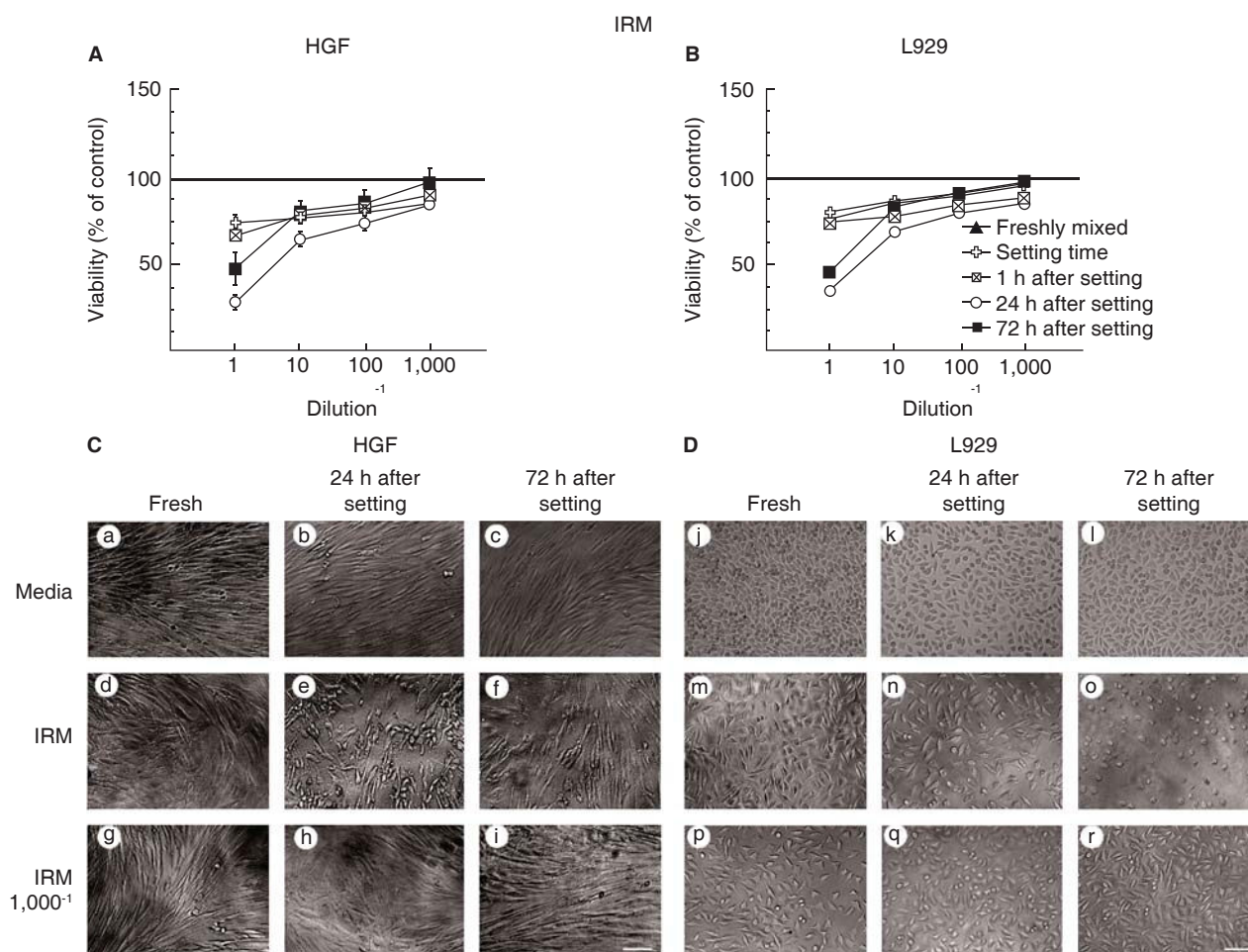


Figure 1. Viability of human gingival fibroblasts (HGF, panel **A**) and L929 mouse cell line fibroblasts (panel **B**) after exposition to IRM[®]. Symbols represent the following: fresh cement (black triangle) setting (white cross) at 1 h (stricken square) at 24 h (white circle) and at 72 h after setting (black square); 100% cell viability represents the state without cement (black line in 100%). Data represent average percentage of cell viability in three independent experiments conducted in triplicate. Microphotographs of HGF (panel **C**) and L-929 (panel **D**) exposed to IRM[®] concentrated conditioning medium and at 1/1,000 diluted IRM[®] conditioning medium. Scale: 100 μ m.

of the ERRM conditioned medium showed higher metabolic activity and greater viability than untreated fibroblasts (Figures 3A and B). In agreement, both human gingival fibroblasts and L929 fibroblast had a typical morphology and remained attached to the culture dish at all times (Figures 3C and 3D, respectively).

DISCUSSION

A yearly estimation of orthodontic procedures mentions that approximately 5.5% of these procedures involve periapical surgery.¹³ In this case, it was necessary to use biocompatible retrofilling materials. Different retrofilling materials exist in the

market, but few comparative reports about their cytotoxicity potential using human gingival fibroblast are available.¹³ For the aforementioned reason, in the present study, we assessed *in vitro* biocompatibility of three retrofilling materials used in periapical surgery: one of them, EndoSequence[®] BM RRM[™] (ERRM, Root Repair Material), recently launched in the Mexican market and two well know materials: IRM[®] and white MTA Angelus[®].

Cytotoxicity was evaluated following the international standard ISO 10993-5 (ISO 10993, 2009), which describes the tests for assessing cytotoxicity *in vitro*. Cytotoxicity was assessed using human gingival fibroblasts (HGF) and L929 mouse fibroblast cell line. Human gingival fibroblasts were

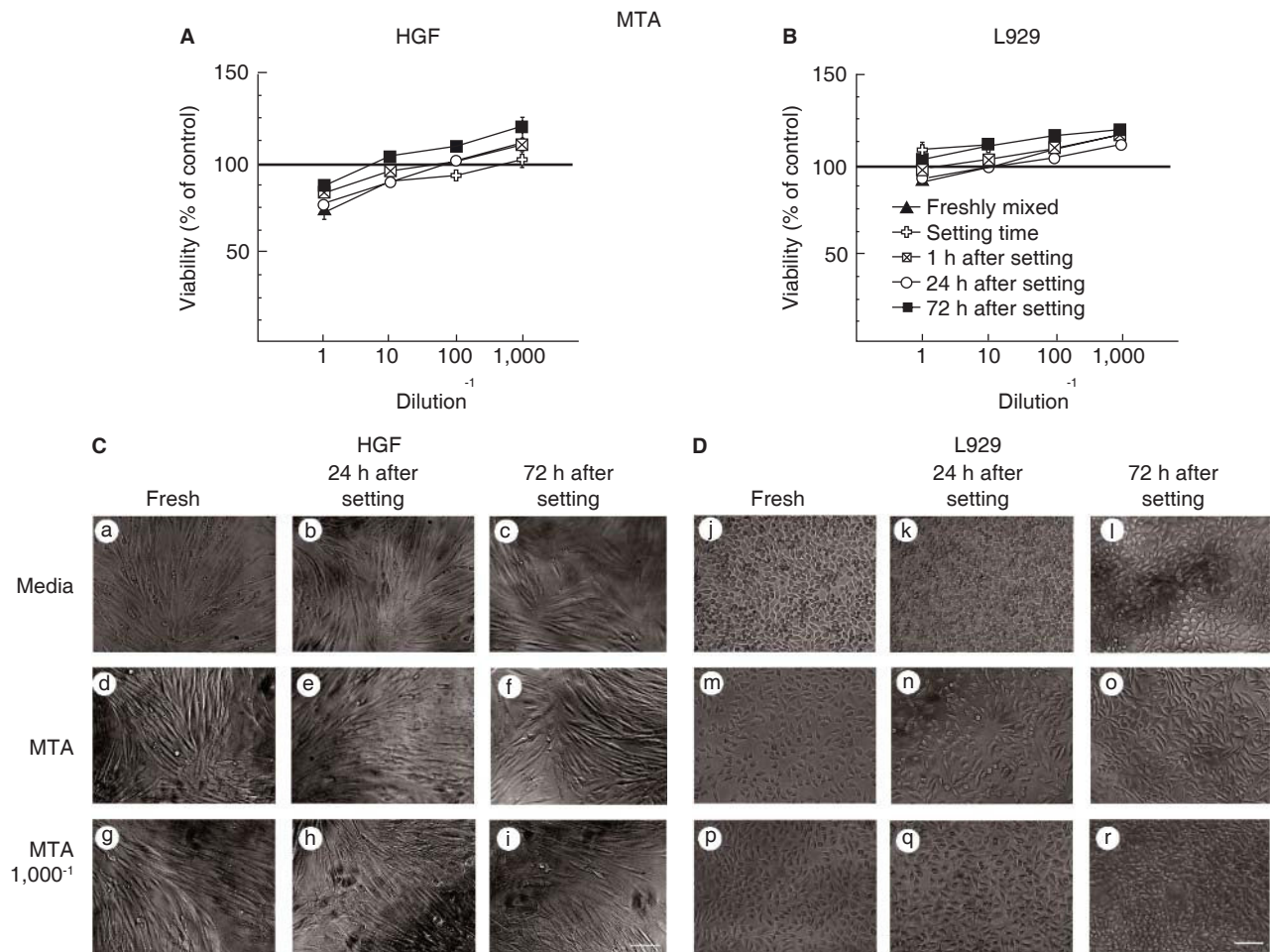


Figure 2. Viability of human gingival fibroblasts (HGF panel A) and mouse fibroblasts of L929 cell line (panel B) after exposition to MTA. Symbols represent the following: fresh cement (black triangle) setting (white cross), after 1 h (stricken square), after 24 h (white circle) and after 72 h (black circle); 100% cell viability represents the state without cement (100% black line). Data represent cell viability percentage average of three independent experiments, conducted in triplicate. HGF microphotographs (panel C) and L929 (panel D) exposed to medium, concentrated MTA conditioned medium and to a 1/1,000 dilution of MTA conditioned medium. Scale 100 μ m.

selected to imitate cell response induced by cements in periapical tissue after endodontic retrofilling. L929 mouse cell line fibroblasts have to their high sensitivity to toxic products.^{22,23} Cytotoxicity was evaluated with the use of two parameters: morphology and cell viability. These parameters were assessed using cements at different circumstances: fresh, set, 12, 24 and 72 hours after setting time.

Most of the studies with cell cultures have focused on the cytotoxicity of materials used in endodontics, especially when the material is set.² When materials are in a fresh state, they release great amounts of chemical sub-products which are toxic to cell cultures. Nevertheless, in virtue of their clinical conditions (*in vivo*), these sub-products dilute in fluids of interstitial

tissues and are eliminated through the vasculature. Therefore, cytotoxic effect of diluting the materials has also being studied.^{2,24}

In the present study we found that IRM[®], a material composed of zinc oxide and eugenol, showed to be the most cytotoxic 24 and 72 hours after set. This fact concurs with other recent studies.^{4,8,25,26} One explanation for this high cytotoxicity might be illustrated by the fact that both eugenol and zinc are cytotoxic.^{18,27} It is thought that released zinc is partly responsible for its prolonged cytotoxic effect.²⁸ Another explanation could be that variations in composition of reinforced materials based on zinc oxide and eugenol could affect their dissolution rate, causing variations in cytotoxicity.^{4,29} Moreover, in IRM[®], eugenol might

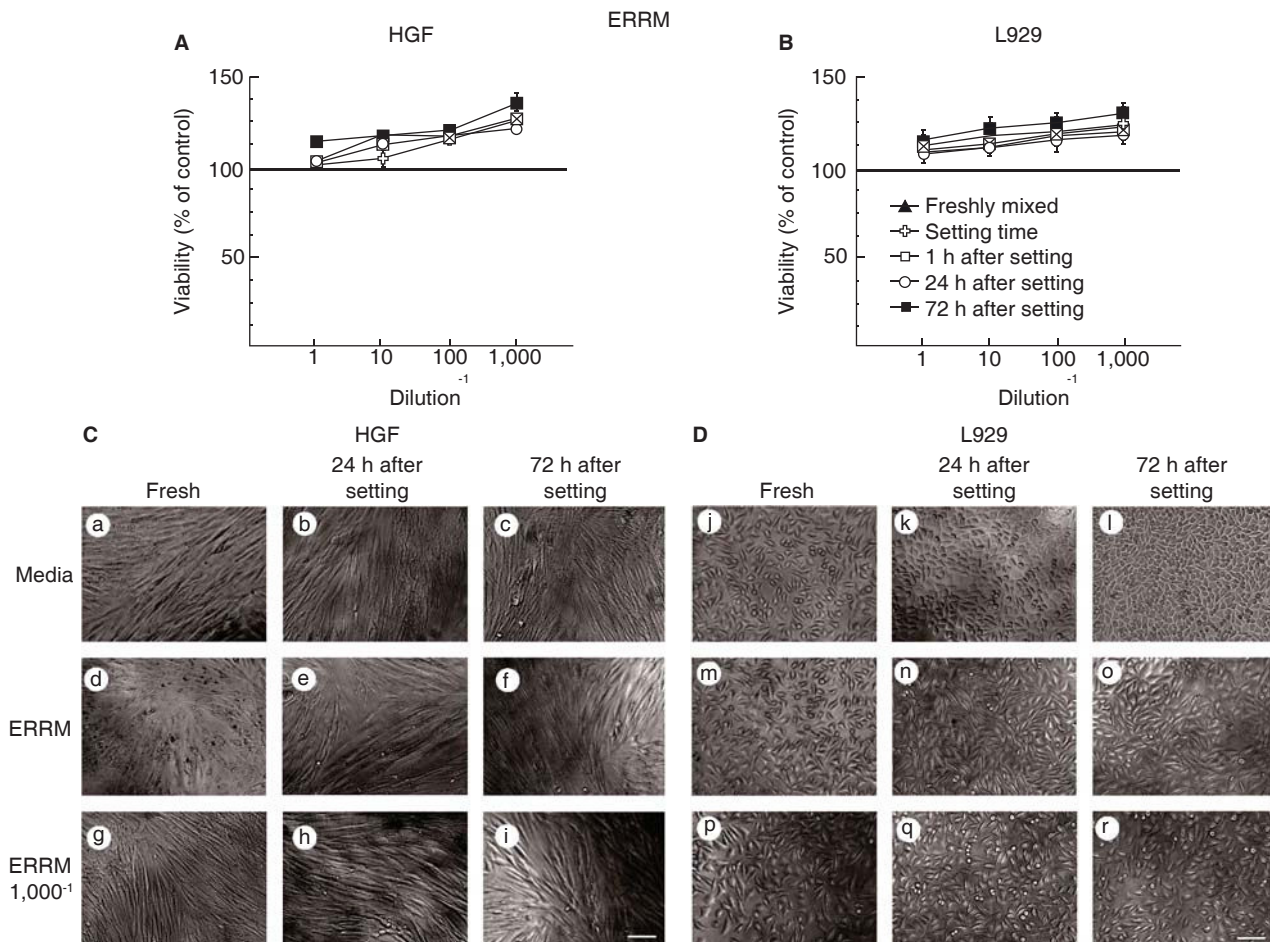


Figure 3. Human gingival fibroblast viability (HGF panel **A**) and fibroblasts of mouse cell line L929 (panel **B**) after exposition to ERRM. Symbols represent the following: fresh cement (black triangle), set (white cross), after 1 h (stricken square) after 24 h (white circle) and 72 h setting (black square); 100% cell viability represents the state without cement (black line 100%). Data represent cell viability percentage average of three independent experiments, conducted in triplicate-HGF microphotographs (panel **C**) and L929 (panel **D**) exposed to medium, to ERRM concentrated conditioned medium and to a 1/1,000 dilution of ERRM conditioned medium. Scale 100 μ m.

exhibit affinity to polymethylmethacrylate, limiting the release of this material and rendering it less cytotoxic.⁴ Another explanation could be that when IRM[®] is used as retrofilling material, it is recommended to use greater proportion of powder than liquid, since this provides great handling advantages, short setting time and decrease of toxicity and solubility.³⁰

Nevertheless, this strong cytotoxicity was not seen if IRM[®] was diluted. In concordance, Trope et al³¹ in a histological study confirmed the favorable response of tissues to IRM[®].⁴ The advantages of IRM[®] are the following: low cost, ease of mixing and manipulation. Several studies have reported positive results when using this material in periradicular surgery.^{3,10,11,31}

Due to the fact that materials used in Endodontics do not possess all these ideal characteristics, MTA

was initially developed as a retrofilling material, and was later used for other procedures such as pulp capping, pulpotomy, apexogenesis, apexification, root perforation reparations, as well as a filling material for root canals. MTA has been recognized as a bioactive, osteoconductor and biocompatible material. Several reviews on MTA's chemical properties, biocompatibility and clinical applications have been published.^{7,13,14} The present study revealed that with mineral trioxide aggregate (MTA) cytotoxicity was almost nil in its fresh state, this concurs with other recent published studies,^{1,6,8,16,32} where MTA has been deemed as a highly biocompatible retrofilling material. Likewise, Kim and Kratchman⁵ mention in their article that MTA is the most biocompatible retrofilling material, and that it can be used with predictable results in endodontic

surgery procedures.³³ In our study, results concurred with others which showed MTA's high biocompatibility in comparison to IRM[®].^{8,17}

The recently launched in the Mexican market retrofilling material ERRM did not show any degree of cytotoxicity indeed, it increased fibroblast metabolic activity; making it the most biocompatible material of this study. Nevertheless ERRM has shown some cytotoxicity, possessing a degree of cytotoxicity similar to MTA's.^{6,16,27}

CONCLUSION

EndoSequence[®] BC RRM[™] root repair material (ERRM) is the most biocompatible retrofilling material and IRM[®] is more cytotoxic than MTA.

Acknowledgments

Research in the authors' laboratory was supported by grant PAPIIT IA202013 (to EUQ) from Dirección General de Asuntos del Personal Académico, Universidad Nacional Autónoma de México, Mexico. To Nancy Mora Pérez, Master in Basic Biomedical Research, for his technical support.

REFERENCES

- Nair U, Ghattas S, Saber M, Natera M, Walker C, Pileggi R. A comparative evaluation of the sealing ability of 2 root-end filling materials: an *in vitro* leakage study using *Enterococcus faecalis*. *Oral Surg Oral Med Oral Pathol Oral Radiol Endod.* 2011; 112 (2): e74-e77.
- Perinpanayagam H. Cellular response to mineral trioxide aggregate root-end filling materials. *J Can Dent Assoc.* 2009; 75 (5): 369-372.
- Chong BS, Pitt-Ford TR, Hudson MB. A prospective clinical study of mineral trioxide aggregate and IRM when used as root-end filling materials in endodontic surgery. *Int Endod J.* 2003; 36 (8): 520-526.
- Chong BS, Pitt-Ford TR. Root-end filling materials: rationale and tissue response. *Endod Topics.* 2005; 11: 114-130.
- Kim S, Kratchman S. Modern endodontic surgery concepts and practice: a review. *J Endod.* 2006; 32 (7): 601-623.
- Damas BA, Wheeler MA, Bringas JS, Hoen MM. Cytotoxicity comparison of mineral trioxide aggregates and EndoSequence bioceramic root repair materials. *J Endod.* 2011; 37 (3): 372-375.
- Parirokh M, Torabinejad M. Mineral trioxide aggregate: a comprehensive literature review--part I: chemical, physical, and antibacterial properties. *J Endod.* 2010; 36 (1): 16-27.
- Torabinejad M, Parirokh M. Mineral trioxide aggregate: a comprehensive literature review--part II: leakage and biocompatibility investigations. *J Endod.* 2010; 36 (2): 190-202.
- Lee BN, Son HJ, Noh HJ, Koh JT, Chang HS, Hwang IN et al. Cytotoxicity of newly developed ortho MTA root-end filling materials. *J Endod.* 2012; 38 (12): 1627-1630.
- Tawil PZ, Trope M, Curran AE, Caplan DJ, Kirakozova A, Duggan DJ et al. Periapical microsurgery: an *in vivo* evaluation of endodontic root-end filling materials. *J Endod.* 2009; 35 (3): 357-362.
- Lindeboom JA, Frenken JW, Kroon FH, van den Akker HP. A comparative prospective randomized clinical study of MTA and IRM as root-end filling materials in single-rooted teeth in endodontic surgery. *Oral Surg Oral Med Oral Pathol Oral Radiol Endod.* 2005; 100 (4): 495-500.
- Islam I, Chng HK, Yap AU. Comparison of the physical and mechanical properties of MTA and portland cement. *J Endod.* 2006; 32 (3): 193-197.
- Roberts HW, Toth JM, Berzins DW, Charlton DG. Mineral trioxide aggregate material use in endodontic treatment: a review of the literature. *Dent Mater.* 2008; 24 (2): 149-164.
- Camilleri J, Pitt-Ford TR. Mineral trioxide aggregate: a review of the constituents and biological properties of the material. *Int Endod J.* 2006; 39 (10): 747-754.
- Obando-Pereda GA, Torres-Chávez KE, Salas-Beltrán H, Hofling JF. Análisis de la composición química, capacidad de sellado apical y propiedades antimicrobianas del MTA y del cemento Portland. *Endodoncia.* 2009; 27 (3): 111-120.
- Alanezi AZ, Jiang J, Safavi KE, Spangberg LS, Zhu Q. Cytotoxicity evaluation of endosequence root repair material. *Oral Surg Oral Med Oral Pathol Oral Radiol Endod.* 2010; 109 (3): e122-e125.
- Yoshimine Y, Ono M, Akamine A. *In vitro* comparison of the biocompatibility of mineral trioxide aggregate, 4 META/MMA-TBB resin, and intermediate restorative material as root-end-filling materials. *J Endod.* 2007; 33 (9): 1066-1069.
- International Organization for Standardization. ISO 10993-5: biological evaluation of medical devices-Part 5: tests for *in vitro* cytotoxicity. ISO 10993-5. 2009.
- Hirschman WR, Wheeler MA, Bringas JS, Hoen MM. Cytotoxicity comparison of three current direct pulp-capping agents with a new bioceramic root repair putty. *J Endod.* 2012; 38 (3): 385-388.
- Torabinejad M, Hong CU, Pitt-Ford TR, Kettering JD. Cytotoxicity of four root end filling materials. *J Endod.* 1995; 21 (10): 489-492.
- Mossman PB, Young LL. Testing for degrees of color blindness. *Occup Health Saf.* 1983; 52 (8): 49-53, 55.
- Silva EJ, Accorsi-Mendonça T, Almeida JF, Ferraz CC, Gomes BP, Zaia AA. Evaluation of cytotoxicity and up-regulation of gelatinases in human fibroblast cells by four root canal sealers. *Int Endod J.* 2012; 45 (1): 49-56.
- Donadio M, Jiang J, He J, Wang YH, Safavi KE, Zhu Q. Cytotoxicity evaluation of Activ GP and Resilon sealers *in vitro*. *Oral Surg Oral Med Oral Pathol Oral Radiol Endod.* 2009; 107 (6): e74-e78.
- Bonson S, Jeansonne BG, Lallier TE. Root-end filling materials alter fibroblast differentiation. *J Dent Res.* 2004; 83 (5): 408-413.
- Lin CP, Chen YJ, Lee YL, Wang JS, Chang MC, Lan WH et al. Effects of root-end filling materials and eugenol on mitochondrial dehydrogenase activity and cytotoxicity to human periodontal ligament fibroblasts. *J Biomed Mater Res B Appl Biomater.* 2004; 71 (2): 429-440.
- Nakayama A, Ogiso B, Tanabe N, Takeichi O, Matsuzaka K, Inoue T. Behaviour of bone marrow osteoblast-like cells on mineral trioxide aggregate: morphology and expression of type I collagen and bone-related protein mRNAs. *Int Endod J.* 2005; 38 (4): 203-210.
- Ma J, Shen Y, Stojicic S, Haapasalo M. Biocompatibility of two novel root repair materials. *J Endod.* 2011; 37 (6): 793-798.
- Balto HA. Attachment and morphological behavior of human periodontal ligament fibroblasts to mineral trioxide aggregate: a scanning electron microscope study. *J Endod.* 2004; 30 (1): 25-29.
- Chong BS, Owadally ID, Pitt-Ford TR, Wilson RF. Cytotoxicity of potential retrograde root-filling materials. *Endod Dent Traumatol.* 1994; 10 (3): 129-133.
- Crooks WG, Anderson RW, Powell BJ, Kimbrough WF. Longitudinal evaluation of the seal of IRM root end fillings. *J Endod.* 1994; 20 (5): 250-252.

31. Trope M, Lost C, Schmitz HJ, Friedman S. Healing of apical periodontitis in dogs after apicoectomy and retrofilling with various filling materials. *Oral Surg Oral Med Oral Pathol Oral Radiol Endod.* 1996; 81 (2): 221-228.
32. Güven EP, Yalvaç ME, Kayahan MB, Sunay H, Şahin F, Bayırlı G. Human tooth germ stem cell response to calcium-silicate based endodontic cements. *J Appl Oral Sci.* 2013; 21 (4): 351-357.
33. Parirokh M, Torabinejad M. Mineral trioxide aggregate: a comprehensive literature review-part III: Clinical applications, drawbacks, and mechanism of action. *J Endod.* 2010; 36 (3): 400-413.

Mailing address:

Eileen Uribe-Querol

E-mail: euquerol@comunidad.unam.mx